

A unique tip design

- SURFLASH™ enables quick and easy confirmation of catheter placement in the vein.
- The catheter is made of polyurethane, which makes it extremely flexible, kink-resistant and durable.
- The radio-opaque stripe on the catheter allows clear x-ray visibility.
- The catheter's puncture resistance level is the same as our TERUMO®'s conventional IV catheter models.


With SURFLASH™, visual confirmation of catheter placement occurs instantaneously due to 'flashback'. Flashback appears between the specially grooved needle and the catheter, providing medical staff with a higher certainty of proper catheter placement.

The specially grooved needle also improves the ability to visualize flashback in patients where venipuncture is more difficult. And because SURFLASH™ is made of

polyurethane, which is firm at room temperature and softer at body temperatures, patient comfort is enhanced, without compromising proper technique.



New flashback feature

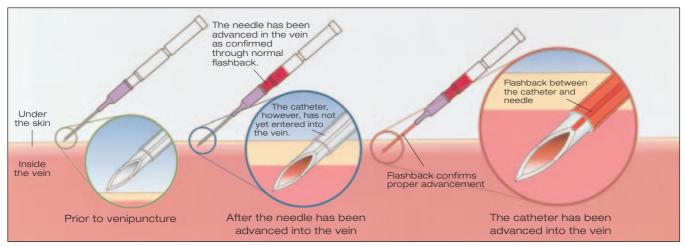

NOTE: The needle insertion angle in the illustrations differs from actual insertion angle.

When using SURFLASH™ for venipuncture, follow the same procedures and take the same safety precautions as when using our TERUMO®'s conventional IV catheter models

References	EN/ISO colour code	External catheter diameter		Internal catheter	Catheter	Flow rate
		Gauge	mm	diameter mm	length mm	ml/min
SR*FF2419	Yellow	24	0.67	0.47	19	15
SR*FF2225	Blue	22	0.85	0.60	25	38
SR*FF2032	Pink	20	1.10	0.80	32	67
SR*FF2051	Pink	20	1.10	0.80	51	59
SR*FF1832	Green	18	1.30	0.95	32	99
SR*FF1851	Green	18	1.30	0.95	51	88
SR*FF1651	Gray	16	1.70	1.30	51	184
SR*FF1464	Orange	14	2.17	1.73	64	282

1 shelf-pack: 50 catheters





Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, Terumo Europe provides you with the clinical systems, means and solutions to keep life flowing.



We keep life flowing

ERUM

LABORATOIRES TERUMO FRANCE S.A. Parc Ariane
Bâtiment Uranus Batiment Uranus Rue Hélène Boucher 78284 Guyancourt Cedex FRANCE Tel. +33/1 30.96.13.00 Fax +33/1 30.43.60.85

TERUMO CORPORATION 44-1, 2-chome Hatagaya

Hatagaya Shibuya-ku Tokyo 151-0072 JAPAN

Tel. +81/3 3374.8111 Fax +81/3 3374.8196

TERUMO EUROPE N.V.

B-3001 Leuven BELGIUM Tel. +32/16.38.12.11 Fax +32/16.40.02.49

BENELUX SALES DIV Interleuvenlaan 40 B-3001 Leuven BELGIUM Tel. +32/16.39.25.80 Fax +32/16.39.25.99 THE NETHERLANDS

Tel. +0800/022.03.96 Fax +0800/022.04.14

TERUMO EUROPE N.V. BENELUX SALES DIVISION

Researchpark Zone 2 Haasrode Interleuvenlaan 40

TERUMO DEUTSCHLAND GmbH

Hauptstrasse 87 65760 Eschborn GERMANY Tel. +49/61.96.80.230 Fax +49/61.96.80.23.200

TERUMO EUROPE ESPAÑA SL Avenida de Burgos 16-D Edificio EUROMOR Madrid 28036 ESPAÑA Tel. +34/91 383.90.69 +34/91 383.91.61 Fax +34/91 767.04.60

TERUMO ITALIA SRL

Via Simone Martini 143/145 I-00142 Roma ITALIA Tel. +39/06.519.61.420 Fax +39/06.503.04.07

TERUMO SWEDEN AB

Nya Varvet Byggnad 90 SE-426 71 Västra Frölunda SWEDEN Tel. +46/31.748.58.80 Fax +46/31.748.58.90

TERUMO UK Ltd.

3 Unity Grove Off School Lane Knowsley Industrial Park South Knowsley Merseyside UNITED KINGDOM Tel. +44/151.482.4800 Fax +44/151.482.4811

®: Registered Trademark Published by Terumo Europe N.V.

www.terumo.com www.terumo-europe.com